


Pocket DAQ manual (rev. 17080)

Apr. 15, 2015

Satoru Yamada


Connection between components (COPPER-II)

About port assignment of FTSW, see
<svn+ssh://xxx@bdaq.local.kek.jp/bdaq/svn/firmware/ft2u/trunk/doc/ft2u.pdf>

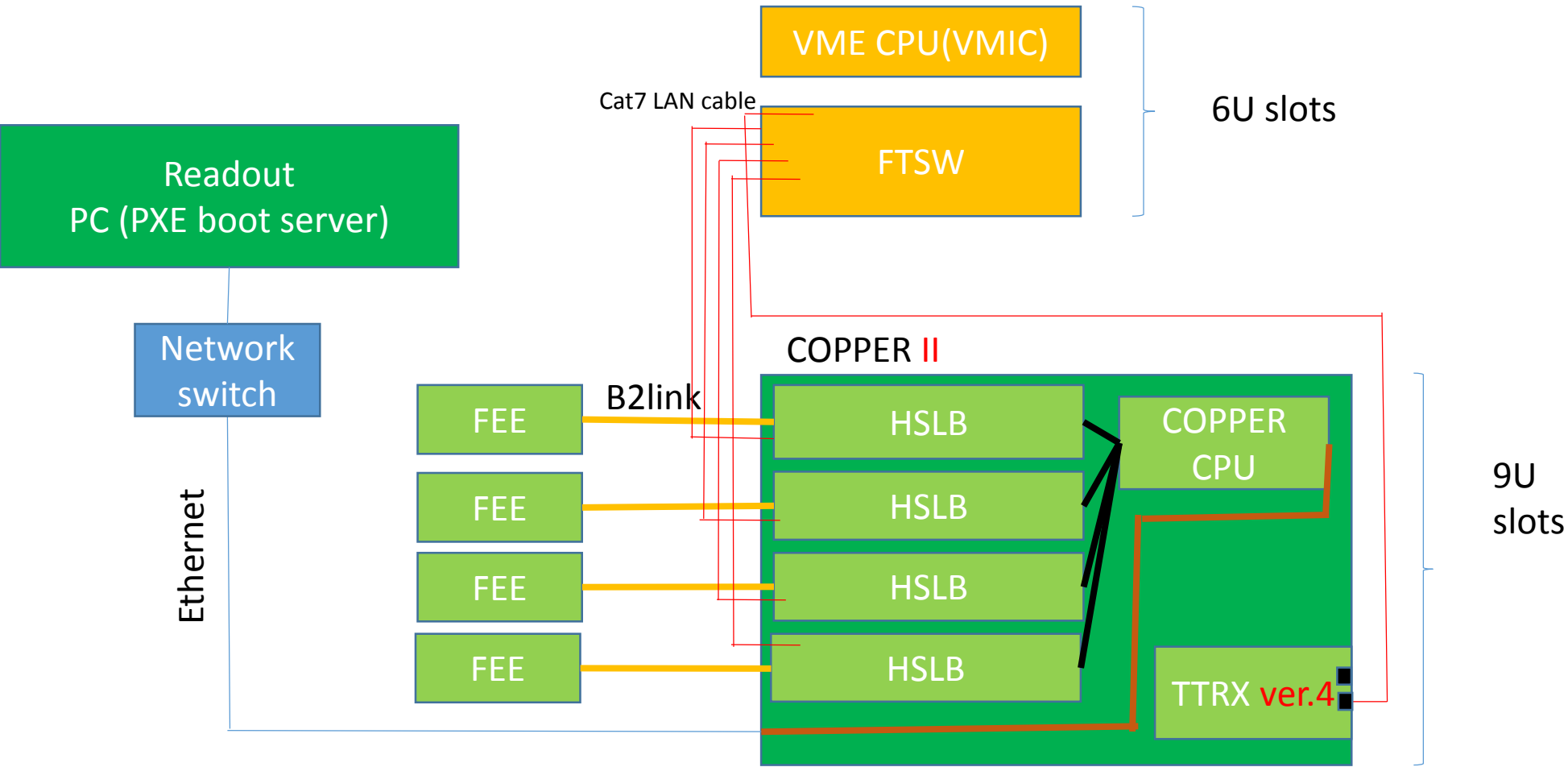

Connection between components (COPPER-III)


About port assignment of FTSW, see
<svn+ssh://xxx@bdaq.local.kek.jp/bdaq/svn/firmware/ft2u/trunk/doc/ft2u.pdf>

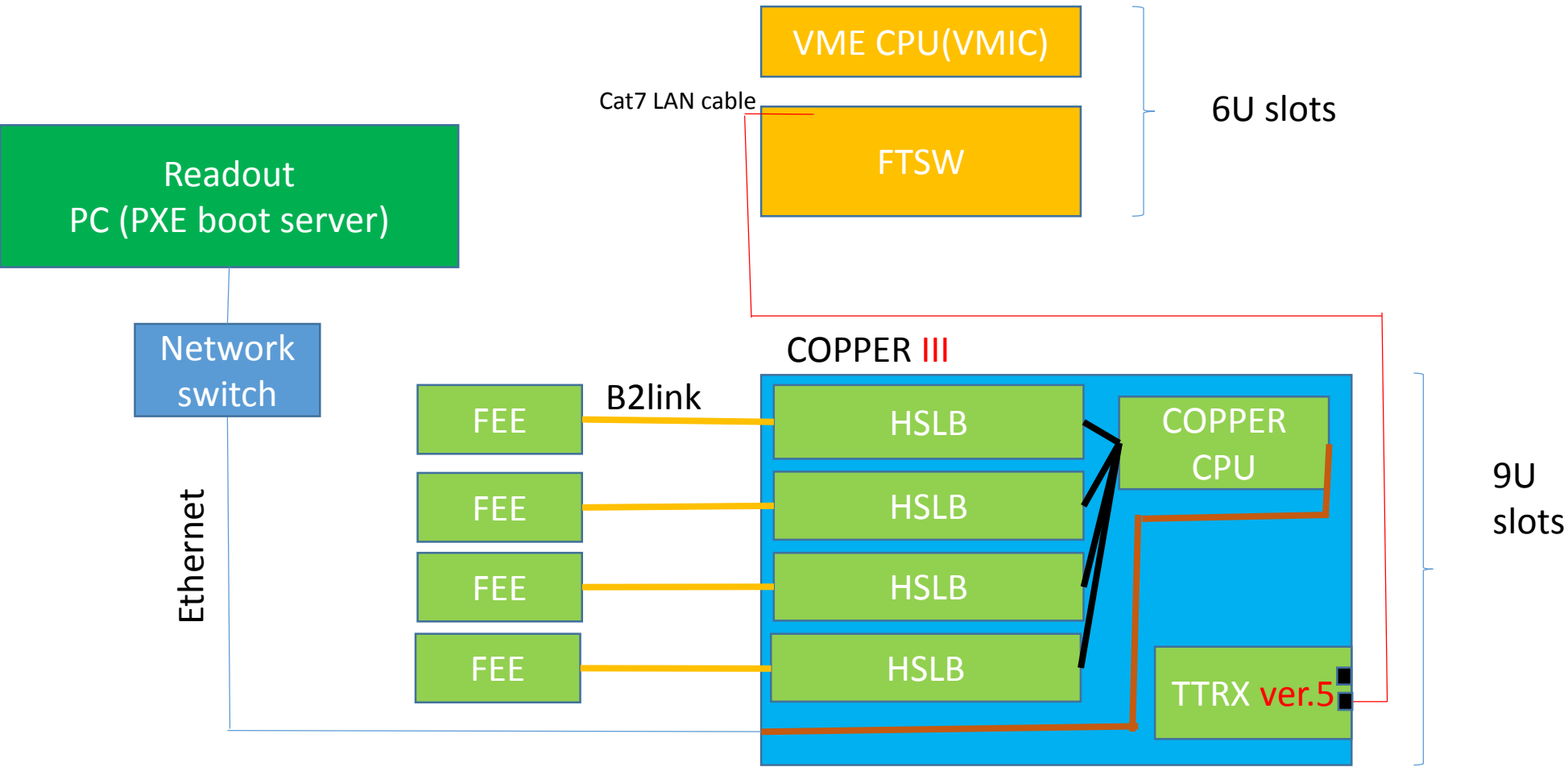


Old Software components

These are stored under <https://belle2.cc.kek.jp/svn/trunk/software/daq>


COPPER CPU

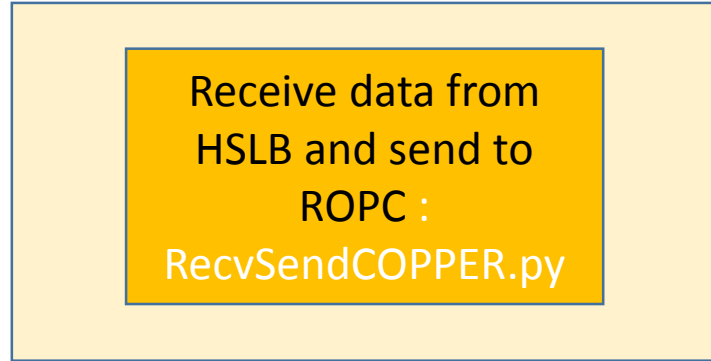

COPPER CPU

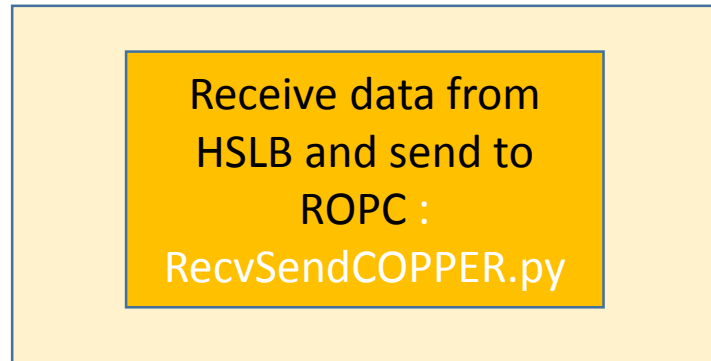

Readout PC

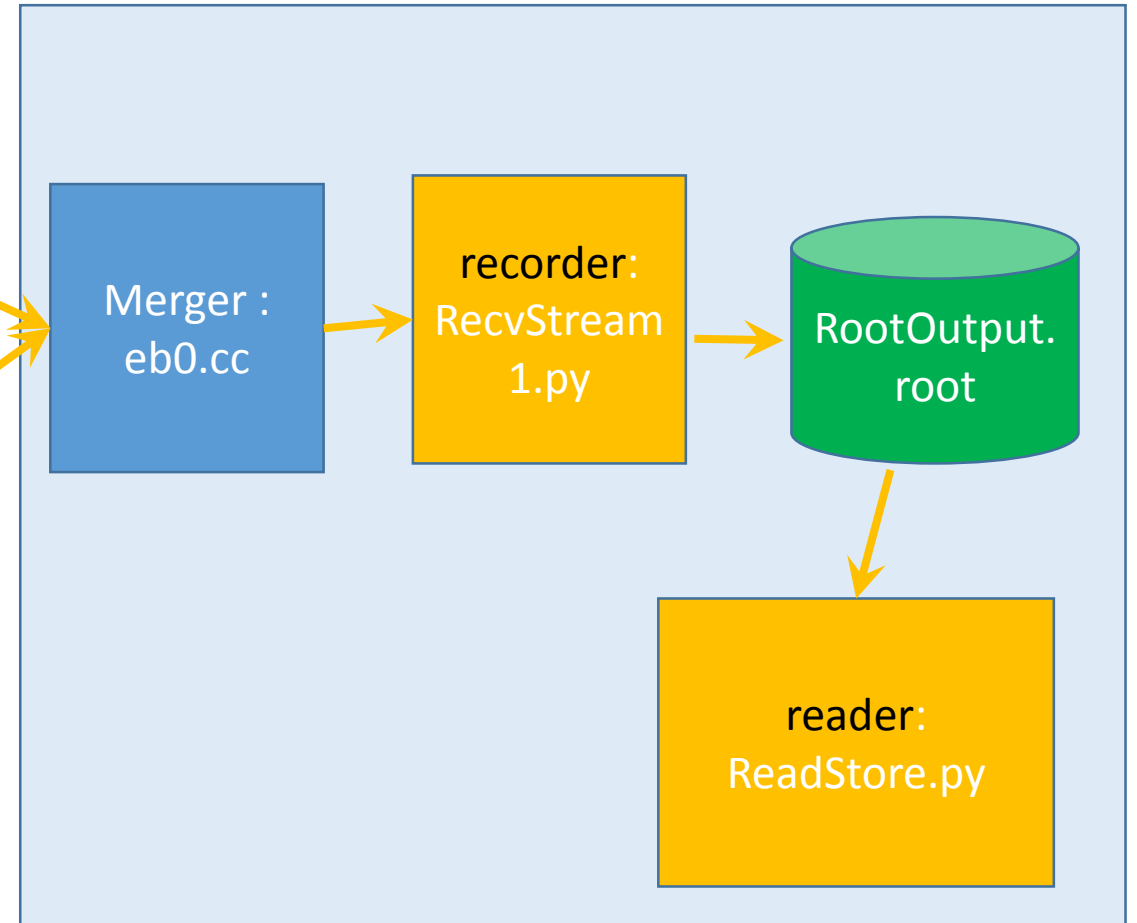


Orange color components are basf2 components.


Old Software components (from rev. 13542 to rev. 15005)

These are stored under <https://belle2.cc.kek.jp/svn/trunk/software/daq>

COPPER CPU

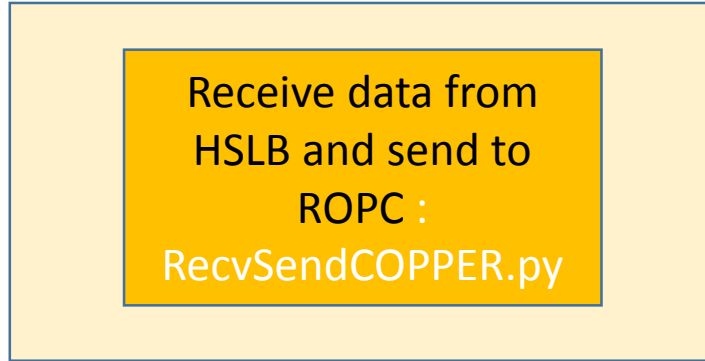


COPPER CPU

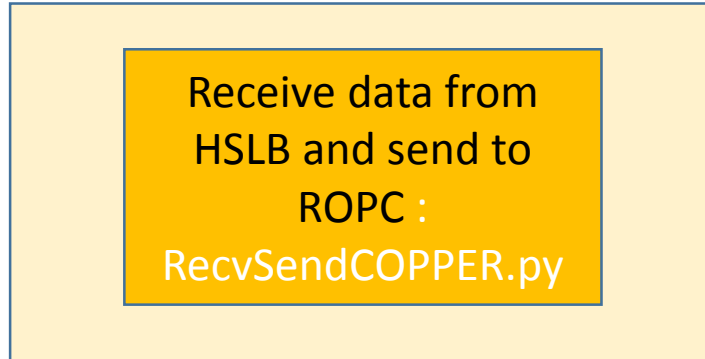

Orange color components are basf2 components.


Readout PC

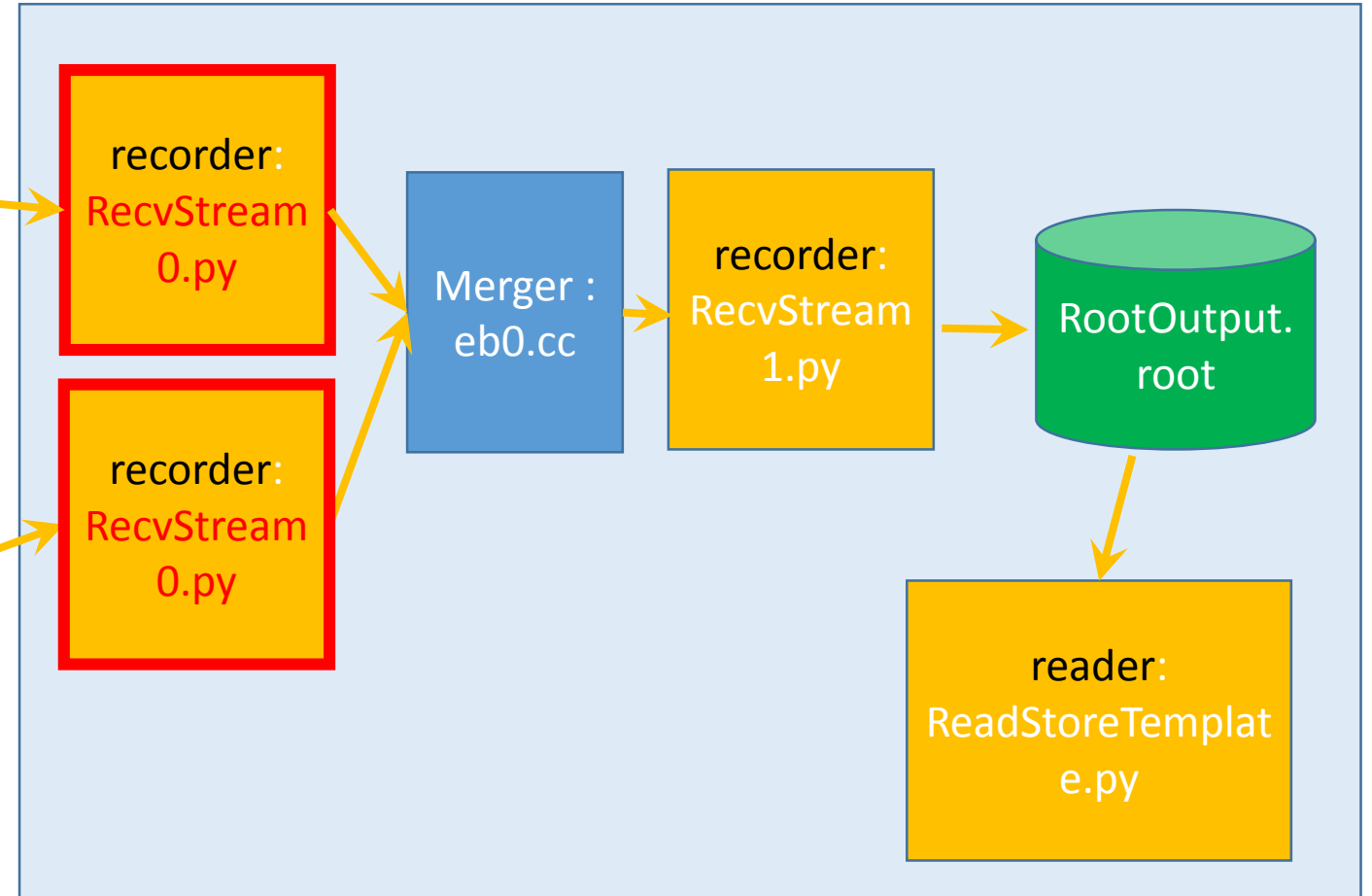


New Software components (from rev. 15006)

These are stored under <https://belle2.cc.kek.jp/svn/trunk/software/daq>


COPPER CPU

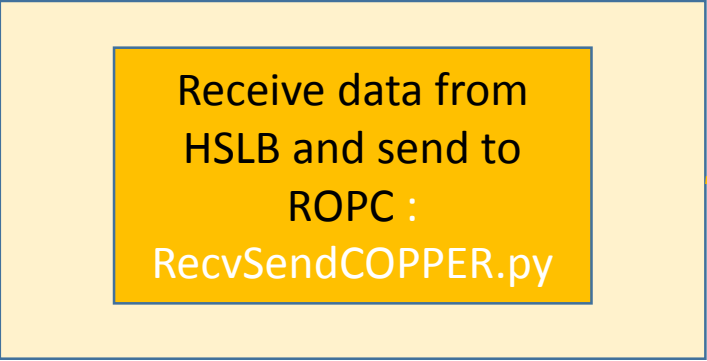

COPPER CPU

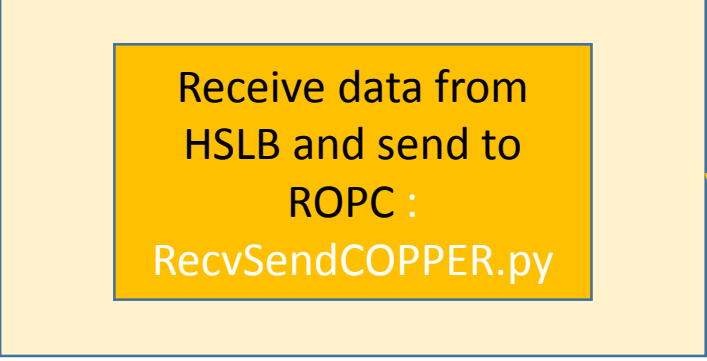

Readout PC

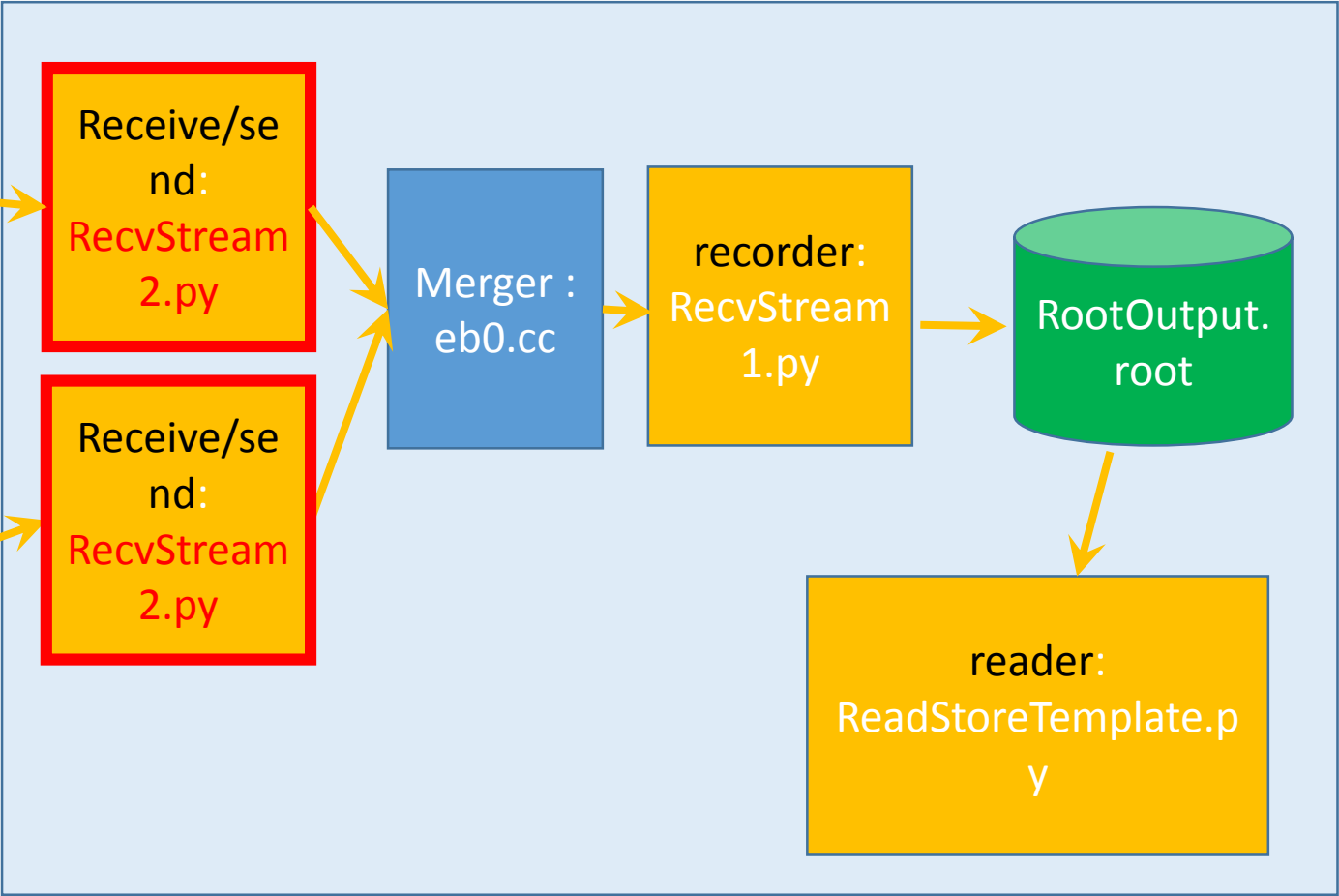

Orange color components are basf2 components.

Pocket DAQ w/o slow
controller and GUI

0. Before using pocket DAQ

- Setup a PXE boot server for COPPER CPU and install driver for COPPER etc
 - See <https://belle2.cc.kek.jp/~twiki/bin/view/Detector/DAQ/PocketDAQ>
- Install basf2 on both COPPER CPU and Readout PC
 - [See https://belle2.cc.kek.jp/~twiki/bin/view/Software/SoftwareInstallation](https://belle2.cc.kek.jp/~twiki/bin/view/Software/SoftwareInstallation) .
“addall” does not install the daq package. You need to do “addpkg daq”.
 - Check daq/Sconscripts and check if `env['CONTINUE'] = False` is commented out.
 - Compile with scon.
- Compile eventbuilder
 - `cd ${BELLE2_LOCAL_DIR}/daq/eventbuilder/evb0/ ; gmake eb0`

Notice from Konno-san about compiling daq package

IMPORTANT!

- Please run the following command to update nsm2 programs in your external directory

```
$ source ${BELLE2_LOCAL_DIR}/daq/slc/extra/nsm2/export.sh
```

Otherwise you may get error messages about compilation of slow control programs ;

daq/slc/nsm/src/NSMData.cc: In member function 'void Belle2::NSMData::parse(const char*, bool)':*

daq/slc/nsm/src/NSMData.cc:299:68: error: too many arguments to function 'NSMparse nsmlib_parsefile(const char*, int, const char*, char*)'*

In file included from include/daq/slc/nsm/NSMData.h:20:0,

from daq/slc/nsm/src/NSMData.cc:1:

Install GLEW for running basf2 program

After external ver.0.5.0?, basf2 requires GLEW(The OpenGL Extension Wrangler Library).

The following procedure worked fine.

1, install the yum epel repository.

```
ropc01$ wget http://ftp-srv2.kddilabs.jp/Linux/distributions/fedora/epel/5/i386/epel-release-5-4.noarch.rpm
```

```
ropc01$ ssh cpr** -lroot
```

```
cpr**# rpm -ivh epel-release-5-4.noarch.rpm
```

```
cpr**# ls -lrt /etc/yum.repos.d/
```

```
cpr**# exit
```

2, install glew on COPPER

```
ropc01$ sudo yum install glew --installroot=/tftpboot/copper/root
```

0.5, Set active HSLB slots in TTRX register

From Nakao-san's mail : [b2link_ml:0159] Re: ft2u-0.65, b2tt-0.31, hslb-0.36

To partially mask HSLB, one can modify register 130 of tt4r:

- setting each of bit 0..3 to 0 corresponds to masking HSLB A..D
- default value of register 130 is 0010000f (all HSLB are enabled)
- bit 20..23 are read only and should be 1 as above
- for example, one can mask HSLB-A by

copper_CPU \$ regrx 130 0000000e

Example :

When you use;

Slot A : **copper_CPU \$ regrx 130 00000001**

Slot B : **copper_CPU \$ regrx 130 00000002**

Slot A and B : **copper_CPU \$ regrx 130 00000003**

2015.4.15

For tt5r028.bit firmware, the setting value should be different.

~~Slot A and B : copper_CPU \$ regrx 130 00000003~~

Slot A and B : copper_CPU \$ regrx 130 00002503

Please ask Nakao-san about details.

1, Set parameters in run_start.sh (1)

```
[ROPC] % cd ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts
```

```
[ROPC] % emacs run_start.sh
```

Set arguments of start_copper.sh

Usage : start_copper.sh **<HOSTNAME>** **<COPPER node>** **<FINNESSE bit flag: A=1, B=2, C=4, D=8>** **<Use network shared memory : Yes=1, No=0>** **<node name(used for Network Shared Memory)>**

e.g.)

```
/usr/bin/xterm -e ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/start_copper.sh cpr006 1 1 0 0&
```

```
/usr/bin/xterm -e ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/start_copper.sh cpr007 2 3 0 0&
```

<COPPER node ID> will be attached to RawCOPPER header.

1, Set parameters in run_start.sh (2)

NOTICE : Environment variables for basf2

In `${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/copper.sh`, a line, “source ~/.bash_profile”, is for setting up basf2 environment. You need to add basf2 setting commands in your `.bash_profile` (or other script file).

Please see "Setup of Software Tools" at <https://belle2.cc.kek.jp/~twiki/bin/view/Software/SoftwareInstallation> for details.

```
[ ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/copper.sh ]
```

```
#
```

```
# setup basf2 environment (See "Setup of Software Tools" at
```

```
https://belle2.cc.kek.jp/~twiki/bin/view/Computing/SoftwareInstallation
```

```
#
```

```
source ~/.bash_profile
```

1, Set parameters in run_start.sh (3)


Set arguments of RecvStream2.sh

Usage : RecvStream2.py <COPPER hostname> <Use NSM(Network Shared Memory)? yes=1/no=0> <port # to accept connection from eb0> <NSM nodename>

Example :

/usr/bin/xterm -e "\${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/RecvStream2.sh cpr001 0 34001 hogehoge1000;" &
/usr/bin/xterm -e "\${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/RecvStream2.sh cpr007 0 34007 hogehoge1000;" &

Readout PC

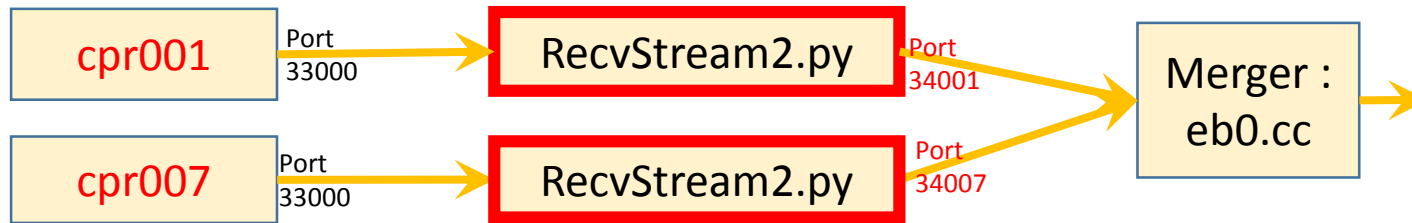

1, Set parameters in run_start.sh (4)

Specify whether communication with GUI via shared memory will be used or not :

1, Edit \${BELLE2_LOCAL_DIR}/daq/rawdata/examples/RecvStream1.py

If you use GUI, set 'UseShmFlag' = 1

```
receiver.param('UseShmFlag', 1)
```

If you don't use GUI set 'UseShmFlag' = 0

```
receiver.param('UseShmFlag', 0)
```

2, Edit \${BELLE2_LOCAL_DIR}/daq/rawdata/examples/RecvStream2.py

If you use GUI, set 'UseShmFlag' = 1

```
receiver.param('UseShmFlag', 1)
```

If you don't use GUI set 'UseShmFlag' = 0

```
receiver.param('UseShmFlag', 0)
```

2, Setting of merger (eb0) (1)

1, compile

```
$ cd ${BELLE2_LOCAL_DIR}/daq/eventbuilder/evb0
```

```
$ make
```


2, edit eb0.sh

`cd /home/usr/b2daq/eb` # change to `cd "your local BelleII directroy"/daq/eventbuilder/evb0/eb-xinetd`
(e.g. `cd /home/usr/yamada/basf2/release/daq/eventbuilder/evb0/eb-xinetd`)

If you want to read data from COPPERs with hostnames of cpr010, cpr011 and cpr012 and corresponding port # of RecvStream2.py are 34010, 34011 and 34012, respectively, use the following ;

~~`./eb0 cpr010:33000 cpr011:33000 cpr012:33000`~~

`./eb0 localhost:34010 localhost:34011 localhost:34012`

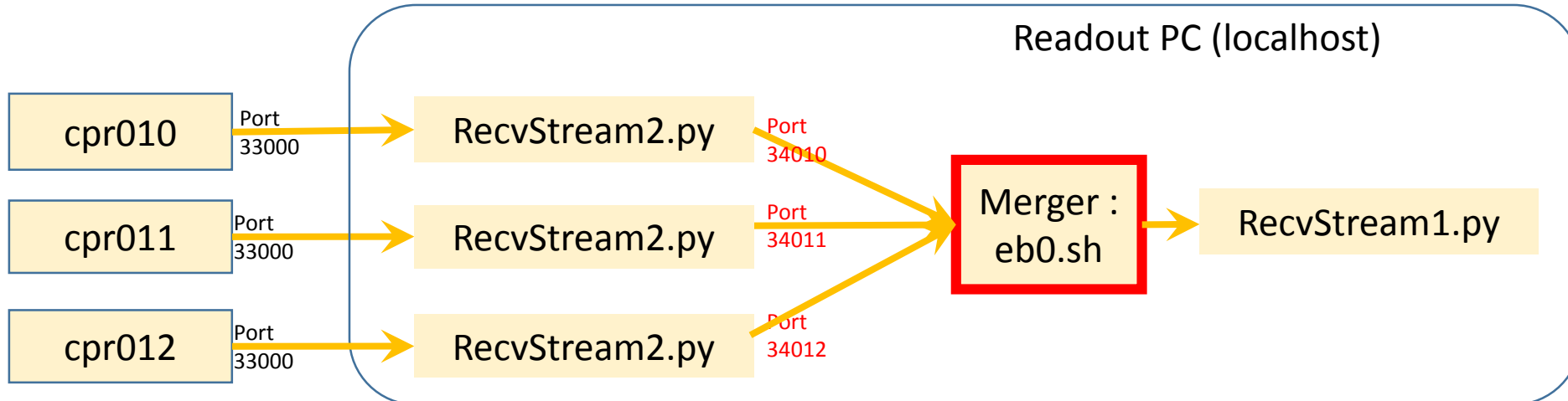

2, Setting of merger (eb0) (2)

3, edit eb0's part in \${BELLE2_LOCAL_DIR}/daq/eventbuilder/evb0/eb-xinted

Modify **user** and **server**.

service eb0

```
{
    bind = 127.0.0.1
    socket_type = stream
    type = UNLISTED
    port = 5101
    wait = no
    user = g0cdc # Change to your user name (e.g. yamagata)
    server = /home/usr/b2daq/eb/eb0.sh
    # change to "your Belle2 local directory"/daq/eventbuilder/evb0/eb0.sh
    # (e.g. /home/usr/yamagata/basf2/release/daq/eventbuilder/evb0/eb0.sh)
}
```

You need to change **the other services in eb-xinetd** (eb1tx, eb1rx, eb2tx, eb2rx) **in the same way.**

2, Setting of merger (eb0) (3)

4, Set link to eb-xinetd

```
# cd /etc/xinetd.d/  
# ln -s ${BELLE2_LOCAL_DIR}/daq/eventbuilder/evb0/eb-  
xinetd eb-xinetd  
# /sbin/service xinetd restart
```

The eb0 daemon will be automatically invoked when a basf2 program on a readout PC connects to a port specified in eb-xinted.

3, How to start DAQ

```
[ROPC] % cd ${BELLE2_LOCAL_DIR}/release/daq/copper/daq_scripts
```

```
[ROPC] % ./run_start.sh
```

It is convenient to login COPPER w/o typing password.

For example;

```
ropc $ cd ~/.ssh
```

```
ropc $ ssh-keygen -t rsa
```

```
ropc $ cat ./id_dsa.pub >> authorized_keys
```

```
ropc $ chmod 600 authorized_keys
```

4, How to stop DAQ

- No stop button for now (If you use GUI, you have.)
- You need to specify max # of events or time to stop the run on a basf2 python file.

```
ROPC % cd ${BELLE2_LOCAL_DIR}/daq/rawdata/examples  
Edit RecvStream1.py
```

You can set following paramters to stop a run.

```
receiver.param('MaxTime', 300.)
```

```
receiver.param('MaxEventNum', 30.)
```

-1 means infinite.

5. Read an output file and extract FEE buffer

- Output file name
 - `${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/root_output.root`
 - You can change the filename by editing `${BELLE2_LOCAL_DIR}/daq/rawdata/examples/RecvStream1.py`
- Read the root file
 - Please see dataformat and unpacker manuals.
 - Data format
 - https://belle2.cc.kek.jp/svn/trunk/software/daq/copper/doc/SetupPocketDAQ_4p1_RawCOPPERDataHandling.pdf
 - Unpacker
 - https://belle2.cc.kek.jp/svn/trunk/software/daq/copper/doc/SetupPocketDAQ_4p5_RawDataUnpackerPacker.pdf

Revision history(1)

- 2013/10/16
 - Change explanation about option of run_start.sh script
- 2014/7/15 (rev. 11619)
 - To specify FINESSE bit field, you need to use decimal value (not hexadecimal).
 - Add an argument of RecvSenDCOPPER.py
 - Update eb0 related issues.
 - About eb0.sh
 - About eb0-xinted
- 2014/7/16(rev. 11654)
 - eb0.sh and eb-xinetd need more modification
- 2014/7/19
 - Add a comment about “addpkg daq”
- 2014/8/8
 - Modify usage of start_run.sh. (nodename is only used for NSM)
- 2014/10/29
 - Change a DAQ scheme so that RecvStream0.py is used on a readout PC
- 2014/11/10
 - Add figures to indicate connection between COPPER and FTSW
 - Add procedure to ssh-login w/o typing password
- 2014/11/14
 - Add “Notice from Konno-san about compiling daq” page
- 2015/1/22
 - Use RecvStream2.py(=DesSerPrePC.cc) instead of RecvStream0.py(=DeSerializerPrePC.cc+Serializer.cc)

Revision history(2)

- 2015 4/13
 - Add “0.5, Set active HSLB slots in TTRX register”

Appendix

A-1, yum setting for COPPER

Sometime yum fails when you try to install software to COPPER directory;

```
readoutPC # yum install **** --installroot=/tftpboot/copper/root
```

In that case, you need to change yum setting for COPPER host:

1, Copy yum repository files.

```
readoutPC # cp /etc/yum.repo.d/* /tftpboot/copper/root/etc/yum.repo.d/*
```

2, Set a proxy server in /tftpboot/copper/root/etc/yum.conf **same as** readoutPC:/etc/yum.conf

A-2, Socket parameters on COPPER, readout PC and etc..

1, Enable "Jumbo Frame" in network switch's setting, if necessary (e.g. HP1810-24G)

If the switch does not support Jumbo frame, setting MTU a large value may cause a problem.

2, Change socket parameters

```
# /sbin/ifconfig eth? mtu 9000
```

```
# /sbin/sysctl -w net.ipv4.tcp_rmem="8388608 8388608 8388608"
```

```
# /sbin/sysctl -w net.ipv4.tcp_wmem="8388608 8388608 8388608"
```

```
# /sbin/sysctl -w net.ipv4.tcp_mem="9000000 9000000 9000000"
```

```
# /sbin/sysctl -w net.core.rmem_default=8388608
```

```
# /sbin/sysctl -w net.core.wmem_default=8388608
```

```
# /sbin/sysctl -w net.core.rmem_max=8388608
```

```
# /sbin/sysctl -w net.core.wmem_max=8388608
```

3, edit /etc/rc.local(or /tftpboot/copper/root/etc/rc.local) if you want the parameters automatically set at boot time

A-3, Use ntpd to adjust time on COPPER

1, Edit /tftpboot/copper/root/etc/rc.d/init.d/ntpd

```
# Source networking configuration.
```

```
. /etc/sysconfig/network
```

```
ntpdate 192.168.10.1
```

```
if [ -f /etc/sysconfig/ntpd ];then
```

```
    . /etc/sysconfig/ntpd
```

```
fi
```

2, check ntpd is started at boot time

```
cpr*** $ chkconfig --list | grep ntp
```

A-4, syslogd

Install syslogd :

```
readoutPC # yum install syslogd --installroot=/tftpboot/copper/root
```

Add entry to /tftpboot/copper/root/etc/syslogd.conf if needed:

```
# Log anything (except mail) of level info or higher.
```

```
# Don't log private authentication messages!
```

```
*.info;mail.none;news.none;authpriv.none;cron.none    /var/log/messages
```

B-1, Test bench at Tsukuba B3

