
# Pocket DAQ manual (rev.13542)

Nov. 14, 2014


Satoru Yamada

### Connection between components (COPPER-II)

About port assignment of FTSW, see svn+ssh://xxx@bdaq.local.kek.jp/bdaq/svn/firmware/ft2u/trunk/doc/ft2u.pdf

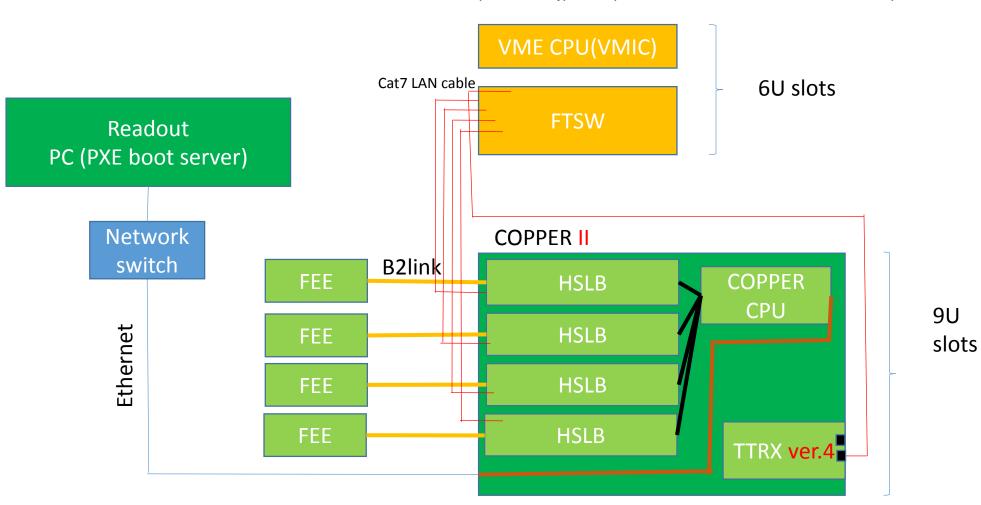


### Connection between components (COPPER-III)

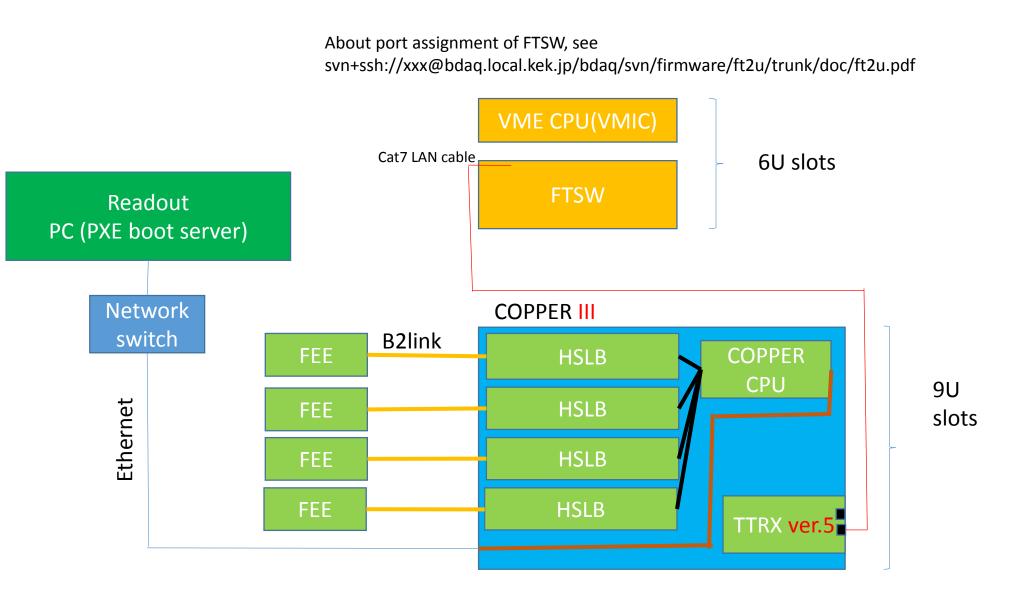

### Old Software components


These are stored under https://belle2.cc.kek.jp/svn/trunk/software/daq

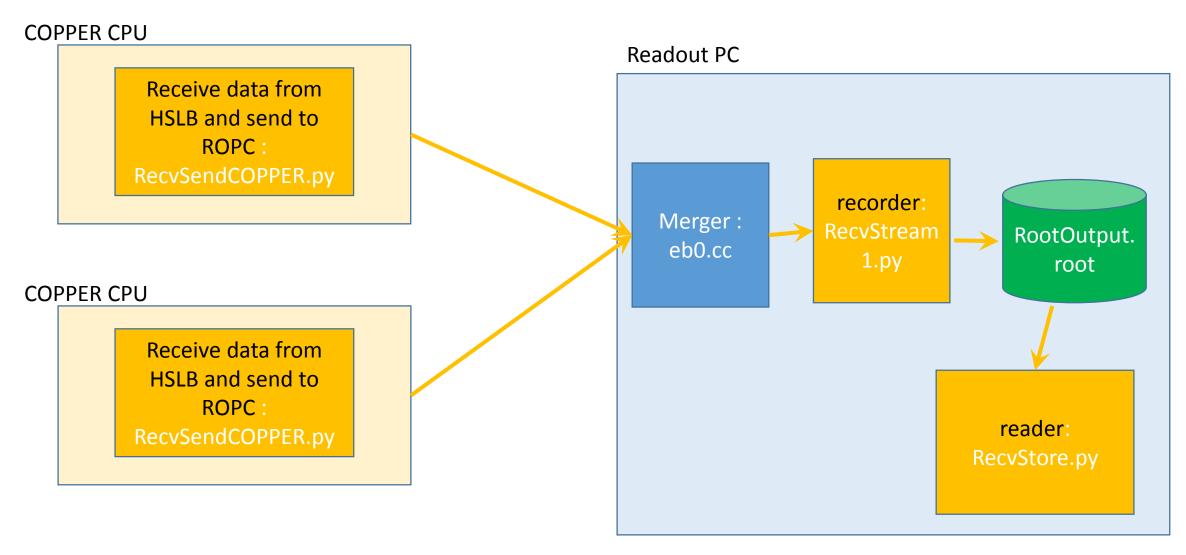

Orange color components are basf2 components.

### New Software components (from rev. 13542)

These are stored under <a href="https://belle2.cc.kek.jp/svn/trunk/software/dag">https://belle2.cc.kek.jp/svn/trunk/software/dag</a>

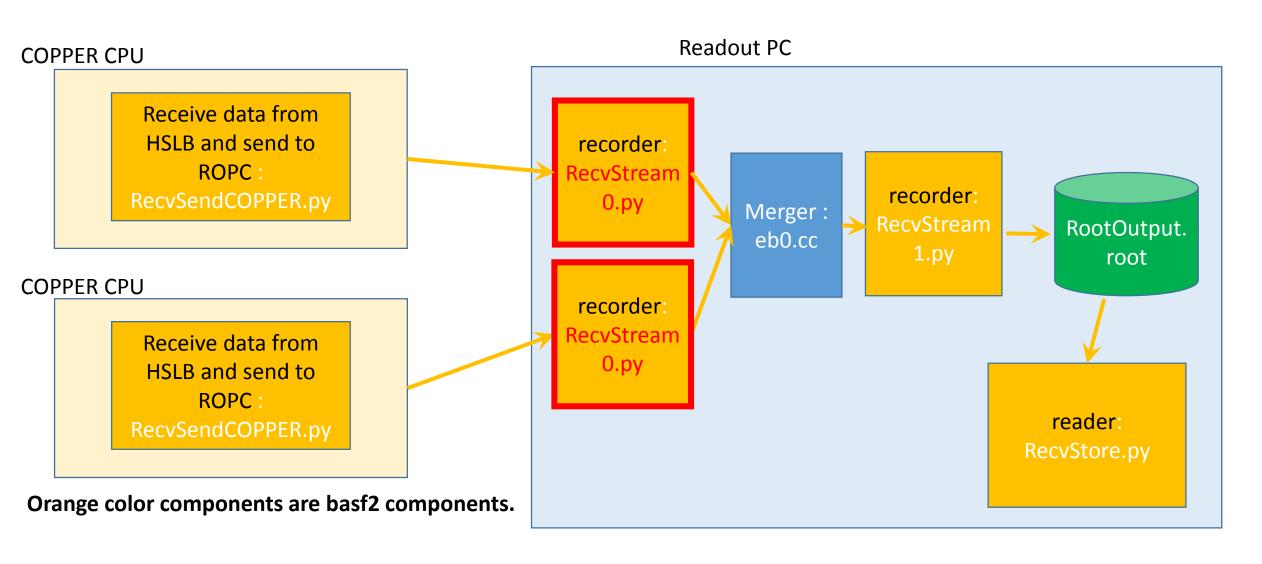

# Pocket DAQ w/o slow controller and GUI

# O. Before using pocket DAQ

- Setup a PXE boot server for COPPER CPU and install driver for COPPER etc
  - See <a href="https://belle2.cc.kek.jp/~twiki/bin/view/Detector/DAQ/PocketDAQ">https://belle2.cc.kek.jp/~twiki/bin/view/Detector/DAQ/PocketDAQ</a>
- Install basf2 on both COPPER CPU and Readout PC
  - <u>See https://belle2.cc.kek.jp/~twiki/bin/view/Software/SoftwareInstallation</u>. "addall" does no install the daq package. You need to do "addpkg daq".
  - Check daq/Sconscripts and check if env['CONTINUE'] = False is commented out.
  - Compile with scons.
- Compile eventbuilder
  - cd \${BELLE2\_LOCAL\_DIR}/daq/eventbuilder/evb0/; gmake eb0

### Notice from Konno-san about compiling daq package

### **IMPORTANT!**

Please run the following command to update nsm2 programs in your external directory

```
$ source ${BELLE2_LOCAL_DIR}/ daq/slc/extra/nsm2/export.sh
```

Otherwise you may get error messages about compilation of slow control programs;

```
daq/slc/nsm/src/NSMData.cc: In member function 'void* Belle2::NSMData::parse(const char*, bool)': daq/slc/nsm/src/NSMData.cc:299:68: error: too many arguments to function 'NSMparse* nsmlib_parsefile(const char*, int, const char*, char*)'
```

In file included from include/daq/slc/nsm/NSMData.h:20:0,

from daq/slc/nsm/src/NSMData.cc:1:

# Install GLEW for running basf2 program

```
After external ver.0.5.0?, basf2 requires GLEW(The OpenGL Extension Wrangler Library).

The following procedure worked fine.

1, install the yum epel repository.

ropc01$ wget http://ftp-srv2.kddilabs.jp/Linux/distributions/fedora/epel/5/i386/epel-release-5-4.noarch.rpm

ropc01$ ssh cpr** -lroot

cpr**# rpm -ivh epel-release-5-4.noarch.rpm

cpr**# ls -lrt /etc/yum.repos.d/

cpr**# exit

2, install glew on COPPER

ropc01$ sudo yum install glew --installroot=/tftpboot/copper/root
```

# 1, Set parameters in run\_start.sh (1)

```
[ROPC] % cd ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts [ROPC] % emacs run_start.sh
```

### Set arguments of start\_copper.sh

```
Usage: start_copper.sh <hOSTNAME> <COPPER node> <FINNESSE bit flag: A=1, B=2, C=4, D=8> <Use network shared memory: Yes=1, No=0> <node name(used for Network Shared Memory)>
e.g.)
/usr/bin/xterm -e ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/start_copper.sh cpr006 1 1 0 0&
/usr/bin/xterm -e ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/start_copper.sh cpr007 2 3 0 0&
```

<COPPER node ID> will be attached to RawCOPPER header.

# 1, Set parameters in run\_start.sh (2)

#### **NOTICE:** Environment variables for basf2

In \${BELLE2\_LOCAL\_DIR}/daq/copper/daq\_scripts/copper.sh, a line, "source ~/.bash\_profile", is for setting up basf2 environment. You need to add basf2 setting commands in your .bash\_profile (or other script file).

Please see "Setup of Software Tools" at https://belle2.cc.kek.jp/~twiki/bin/view/Software/SoftwareInstallation for details.

```
[ ${BELLE2_LOCAL_DIR}/daq/copper/daq_scripts/copper.sh ]
#
# setup basf2 environment (See "Setup of Software Tools" at
https://belle2.cc.kek.jp/~twiki/bin/view/Computing/SoftwareInstallation
#
source ~/.bash_profile
```


# 1, Set parameters in run\_start.sh (3)

### Set arguments of RecvStream0.sh

Usage : RecvStream0.py <COPPER hostname> <Use NSM(Network Shared Memory)? yes=1/no=0> <port # to accept connection from eb0 > <NSM nodename>

Readout PC

Example :
/usr/bin/xterm -e "\${BELLE2\_LOCAL\_DIR}/daq/copper/daq\_scripts/RecvStream0.sh cpr001 0 34001 hogehoge1000;" & /usr/bin/xterm -e "\${BELLE2\_LOCAL\_DIR}/daq/copper/daq\_scripts/RecvStream0.sh cpr007 0 34007 hogehoge1000;" &

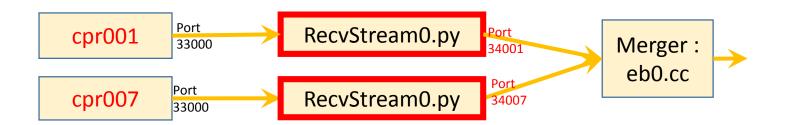

# 1, Set parameters in run\_start.sh (4)

Specify whether communication with GUI via shared memory will be used or not:

1, Edit \${BELLE2\_LOCAL\_DIR}/daq/rawdata/examples/RecvStream1.py

```
If you use GUI, set 'UseShmFlag' = 1
receiver.param('UseShmFlag', 1)
If you don't use GUI set 'UseShmFlag' = 0
receiver.param('UseShmFlag', 0)
```

2, Edit \${BELLE2\_LOCAL\_DIR}/daq/rawdata/examples/RecvStream0.py

```
If you use GUI, set 'UseShmFlag' = 1
receiver.param('UseShmFlag', 1)
If you don't use GUI set 'UseShmFlag' = 0
receiver.param('UseShmFlag', 0)
```

# 2, Setting of merger (eb0) (1)

#### 1, compile


\$ cd \${BELLE2\_LOCAL\_DIR}/daq/eventbuilder/evb0 \$ make

#### 2, edit eb0.sh

cd /home/usr/b2daq/eb # change to cd "your local BelleII directroy"/daq/eventbuilder/evb0/eb-xinetd (e.g. cd /home/usr/yamada/basf2/release/daq/eventbuilder/evb0/eb-xinetd)

If you want to read data from COPPERs with hostnames of cpr010, cpr011 and cpr012 and corresponding port # of RecvStream0.py are 34010, 34011 and 34012, respectively, use the following; ./eb0 cpr010:33000 cpr011:33000 cpr012:33000

./eb0 localhost:34010 localhost:34011 localhost:34012

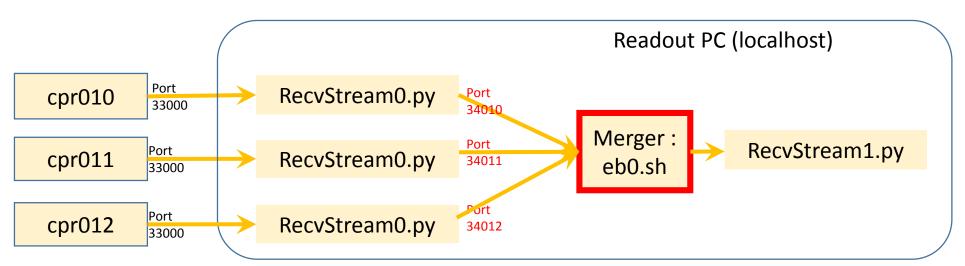

# 2, Setting of merger (eb0) (2)

3, edit eb0's part in \${BELLE2\_LOCAL\_DIR}/daq/eventbuilder/evb0/eb-xinted

```
Modify user and server.

service eb0

{
    bind = 127.0.0.1
    socket_type = stream
    type = UNLISTED
    port = 5101
    wait = no
    user = g0cdc # Change to your user name (e.g. yamagata)
    server = /home/usr/b2daq/eb/eb0.sh
    # change to "your Belle2 local directory"/daq/eventbuilder/evb0/eb0.sh

# (e.g. /home/usr/yamagata/basf2/release/daq/eventbuilder/evb0/eb0.sh)
}
```

You need to change the other services in eb-xinetd (eb1tx, eb1rx, eb2tx, eb2rx ....) in the same way.

# 2, Setting of merger (eb0) (3)

#### 4, Set link to eb-xinetd

```
# cd /etc/xinetd.d/
# In -s ${BELLE2_LOCAL_DIR}/daq/eventbuilder/evb0/eb-
xinetd eb-xinetd
# /sbin/service xinetd restart
```

The eb0 daemon will be automatically invoked when a basf2 program on a readout PC connects to a port specified in ebxinted.

## 3, How to start DAQ

```
[ROPC] % cd ${BELLE2_LOCAL_DIR}/release/daq/copper/daq_scripts [ROPC] % ./run_start.sh
```

```
It is convenient to login COPPER w/o typing password. For example; ropc $ cd ~/.ssh ropc $ ssh-keygen rsa ropc $ cat ./id_dsa.pub >> authorized_keys ropc $ chmod 600 authorized_keys
```

# 4, How to stop DAQ

- No stop button for now (If you use GUI, you have.)
- You need to specify max # of events or time to stop the run on a basf2 python file.

```
ROPC % cd ${BELLE2_LOCAL_DIR}/daq/rawdata/examples Edit RecvStream1.py
```

You can set following paramters to stop a run. receiver.param('MaxTime', 300.) receiver.param('MaxEventNum', 30.)

-1 means inifinite.

## 5. Read an output file and extract FEE buffer

- Output file name
  - \${BELLE2\_LOCAL\_DIR}/daq/copper/daq\_scripts/root\_output.root
  - You can change the filename by editing \${BELLE2\_LOCAL\_DIR}/daq/rawdata/examples/RecvStream1.py
- Read the root file
  - Please see dataformat and unpacker manuals.
  - Data format
    - https://belle2.cc.kek.jp/svn/trunk/software/daq/copper/doc/SetupPocketDAQ\_4p1\_Raw COPPERDataHandling.pdf
  - Unpacker
    - https://belle2.cc.kek.jp/svn/trunk/software/daq/copper/doc/SetupPocketDAQ\_4p5\_Raw DataUnpackerPacker.pdf

### Revision history

- 2013/10/16
  - Change explanation about option of run\_start.sh script
- 2014/7/15 (rev. 11619)
  - To specify FINESSE bit field, you need to use decimal value ( not hexadecimal ).
  - Add an argument of RecvSenDCOPPER.py
  - Update eb0 related issues.
    - About eb0.sh
    - About eb0-xinted
- 2014/7/16(rev. 11654)
  - eb0.sh and eb-xinetd need more modification
- 2014/7/19
  - Add a comment about "addpkg daq"
- 2014/8/8
  - Modify usage of start\_run.sh. (nodename is only used for NSM)
- 2014/10/29
  - Change a DAQ scheme so that RecvStream0.py is used on a readout PC
- 2014/11/10
  - Add figures to indicate connection between COPPER and FTSW
  - Add procedure to ssh-login w/o typing password
- 2014/11/14
  - Add "Notice from Konno-san about compiling daq" page

# Appendix

# A-1, yum setting for COPPER

Sometime yum fails when you try to install software to COPPER directory; readoutPC # yum install \*\*\*\* --installroot=/tftpboot/copper/root

In that case, you need to change yum setting for COPPER host:

- 1, Copy yum repository files. readoutPC # cp /etc/yum.repo.d/\* /tftpboot/copper/root/etc/yum.repo.d/\*
- 2, Set a proxy server in /tftpboot/copper/root/etc/yum.conf same as readoutPC:/etc/yum.conf

### A-2, Socket parameters on COPPER, readout PC and etc..

- 1, Enable "Jumbo Frame" in network switch's setting, if necessary (e.g. HP1810-24G)

  If the switch does not support Jumbo frame, setting MTU a large value may cause a problem.
- 2, Change socket parameters
- #/sbin/ifconfig eth? mtu 9000

```
# /sbin/sysctl -w net.ipv4.tcp_rmem="8388608 8388608 8388608"
# /sbin/sysctl -w net.ipv4.tcp_wmem="8388608 8388608 8388608"
# /sbin/sysctl -w net.ipv4.tcp_mem="9000000 9000000 9000000"
# /sbin/sysctl -w net.core.rmem_default=8388608
# /sbin/sysctl -w net.core.wmem_default=8388608
# /sbin/sysctl -w net.core.rmem_max=8388608
# /sbin/sysctl -w net.core.wmem_max=8388608
```

3, edit /etc/rc.local(or /tftpboot/copper/root/etc/rc.local) if you want the parameters automatically set at boot time

### A-3, Use ntpd to adjust time on COPPER

cpr\*\*\*\* \$ chkconfig –list | grep ntp

### A-4, syslogd

Install syslogd:

readoutPC # yum install syslogd --installroot=/tftpboot/copper/root

### Add entry to /tftpboot/copper/root/etc/syslogd.conf if needed:

# Log anything (except mail) of level info or higher.

# Don't log private authentication messages!

\*.info;mail.none;news.none;authpriv.none;cron.none /var/log/messages

### B-1, Test bench at Tsukuba B3

